

Medtronic

It's just right

Steerant™ Aortic Guidewire

Meet Steerant, the guidewire designed to be just right for your aortic procedures. The guidewire's soft, atraumatic tip gradually transitions to a stiff main body, providing protection for fragile aortic anatomy along with the support needed for device deployment. Steerant guidewire is offered in lengths tailored for EVAR and TEVAR procedures.

Tailored for aortic procedures

- 200 cm wire length is specifically tailored for EVAR procedures
- 260 cm double curved wire to support TEVAR cases
- 260 cm curved wire option available

Protects fragile aortic anatomy

- 15 cm flexible, soft tip provides an atraumatic interface


Balanced control

- Right balance of stiffness and trackability
- Gradual stiffness transition provides balanced control in straightforward and challenging anatomy
- PTFE coating provides smooth device delivery

Right length.
Right stiffness.
Right control.

It's just right.



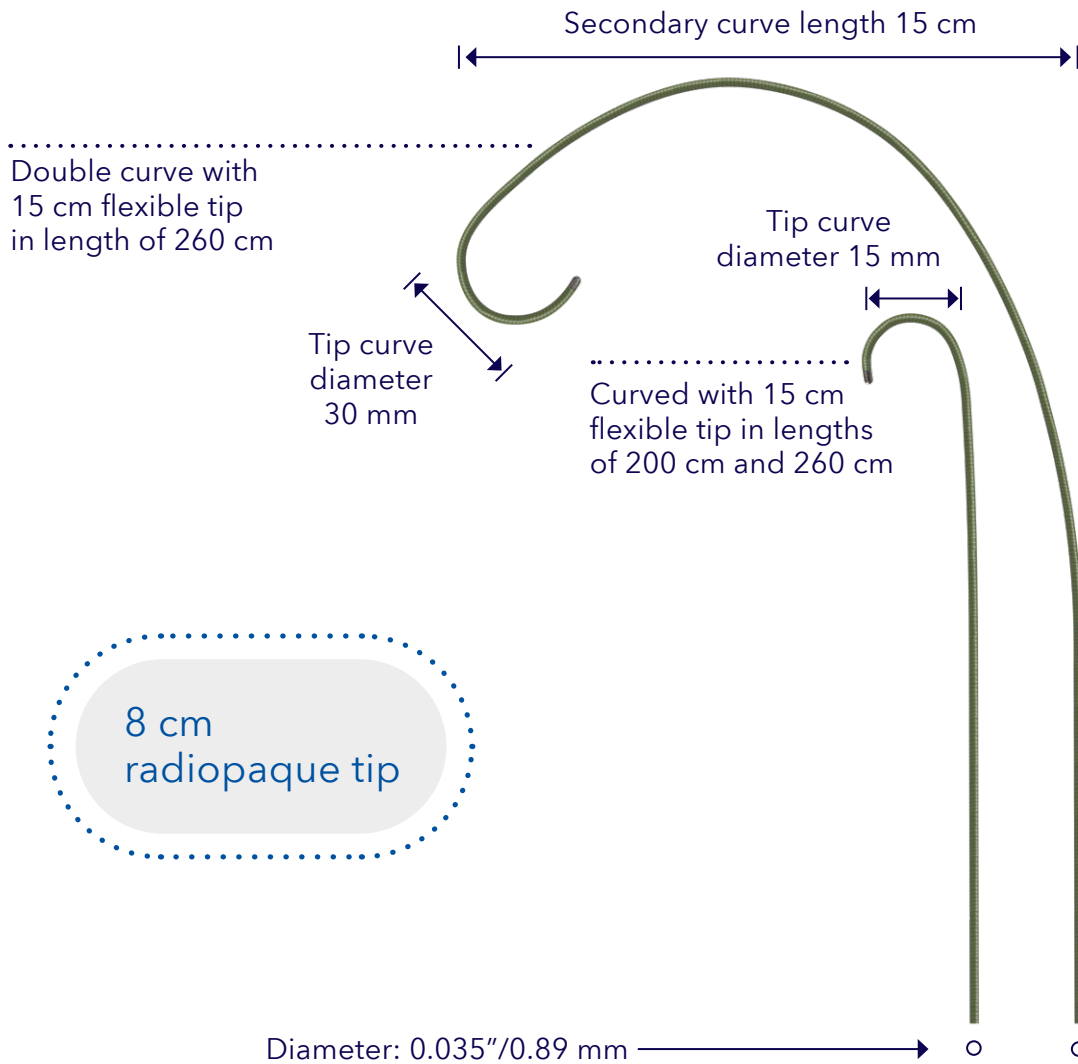

Ordering information

Catalog number	Lengths (cm)	Tip shape	Transition length (cm)	Tip curve diameter (mm)	Secondary curve length (cm)	Wire diameter (in) (mm)	RO length (cm)
AGWSJ200	200	Curved	15	15	N/A	0.035 (0.89)	8
AGWSJ260	260	Curved	15	15	N/A	0.035 (0.89)	8
AGWDJ260	260	Double curve	15	30	15	0.035 (0.89)	8

Important Information: Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device.

Indications for Use: The Medtronic Steerant super stiff guidewire is intended to be used to facilitate catheter placement and exchange during diagnostic or interventional procedures in the aorta. The guidewire is not indicated for coronary or neuro vasculature.

CAUTION: Federal (USA) law restricts these devices to sale by or on the order of a physician.

Medtronic

UC202301879a EN ©2024 Medtronic. Medtronic and the Medtronic logo are trademarks of Medtronic. All other brands are trademarks of a Medtronic company. For distribution in the USA only. 05/2024

[medtronic.com/Steerant](https://www.medtronic.com/Steerant)